228/229 Barrier-Free® Power Examination Chair Setup Instructions

WARNING

Heavy object. Use lifting aids and proper lifting techniques when moving chair to minimize muscle strains or back injury.

Tools:

Battery Powered Drill/Impact Driver
Standard Wrench Set
Standard Deep Well Socket Set
Box Cutter
Furniture Dolly (4 Wheel)


Equipment Alert

When moving the chair, you must adhere to the following instructions.

ALWAYS: a) Transport chair in full-down position with base shroud supports in place.


NEVER: a) Push on base shrouds.
b) Lift chair by stirrups, back section or upholstered top.

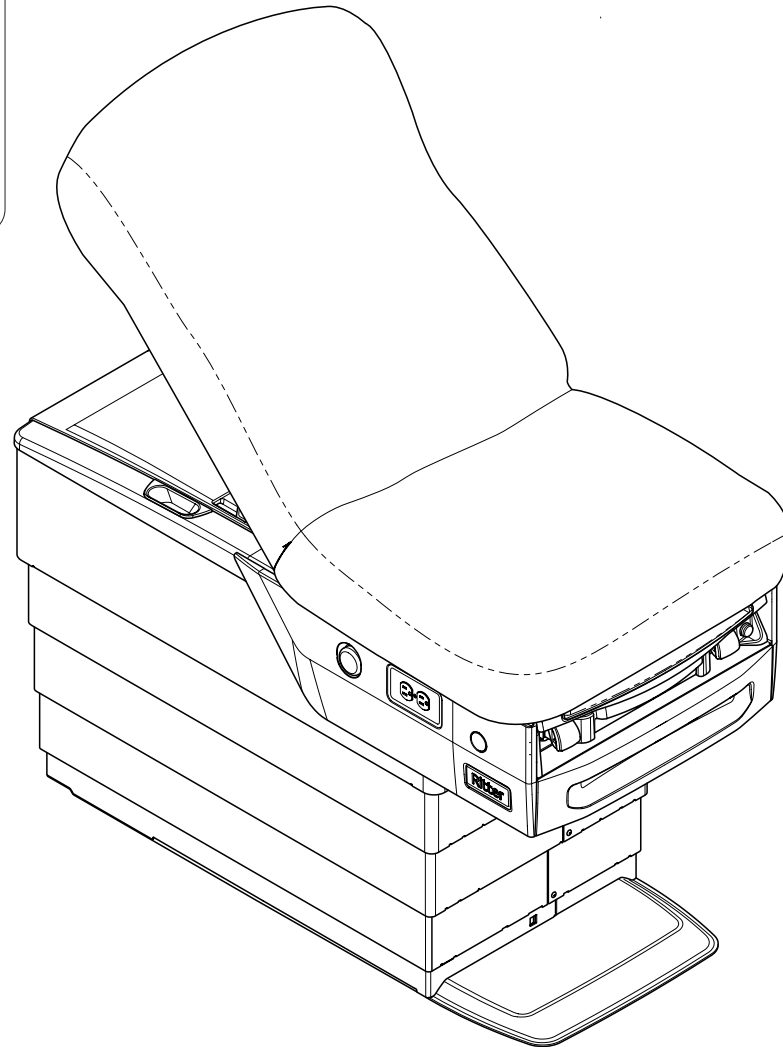
0 0 3 - 1 0 8 8 1 - 9 9


Step 1:

Remove brackets (4) from skid to release chair.

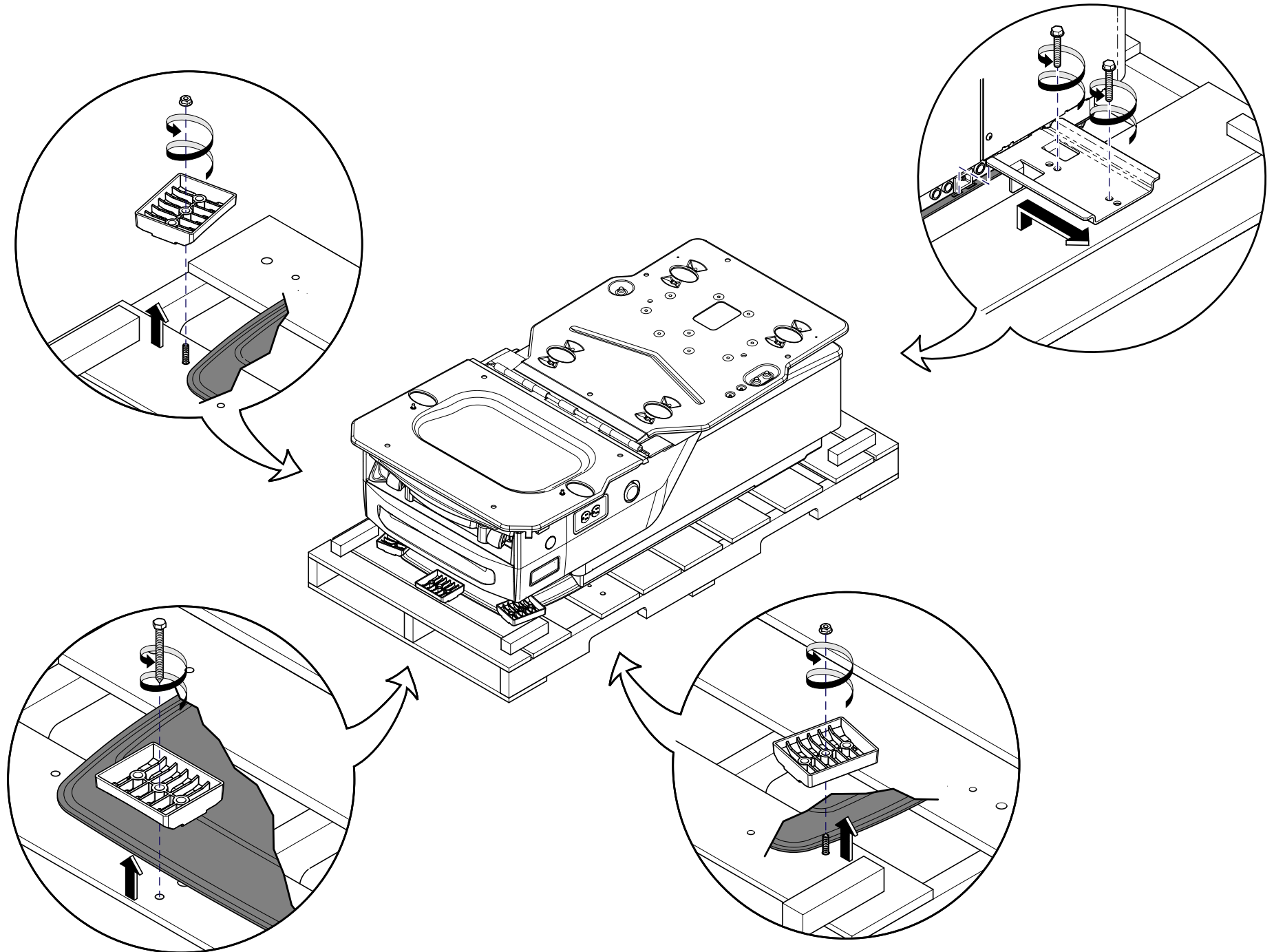
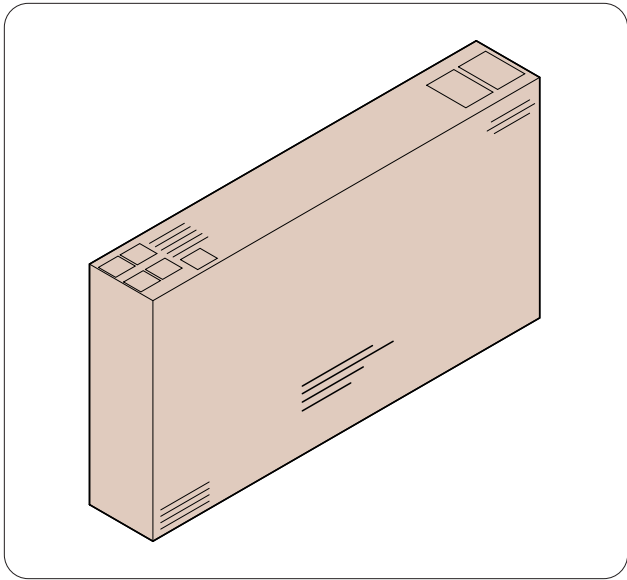
Step 2:
Locate box containing upholstery.

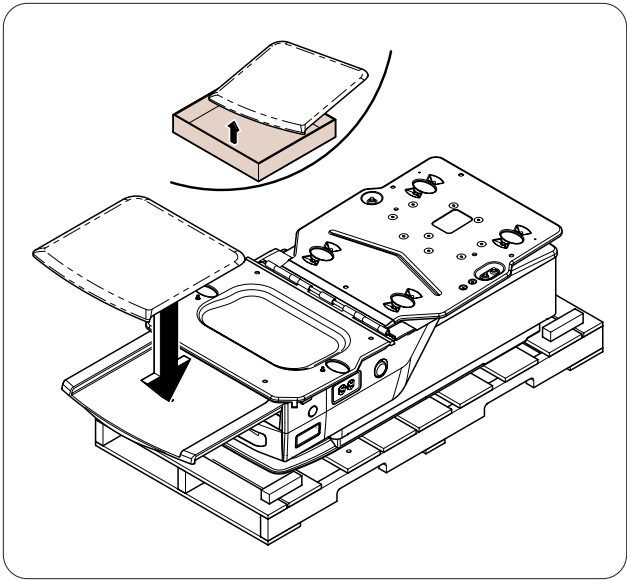
Step 3:
Place footrest upholstery, inserting pin on bottom of upholstery into footrest shelf hole.

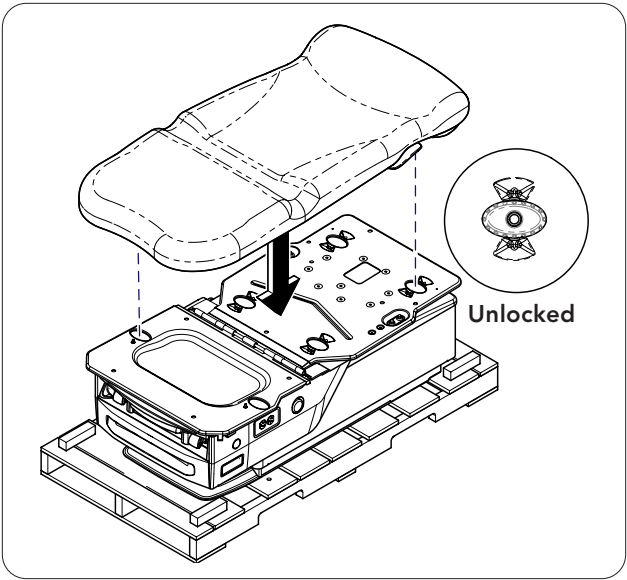
Step 4:
With locking knobs (4) in the unlocked position, place upholstery on chair top aligning locking knobs with holes.

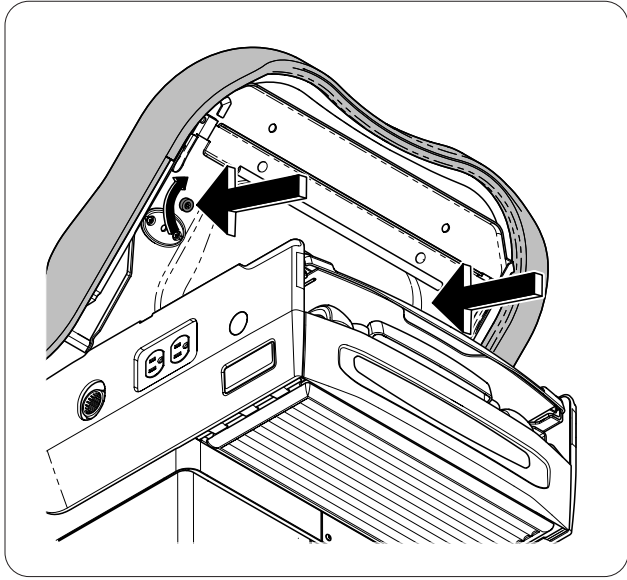
Step 5:
Raise seat pan 90 degrees and secure seat section using attached thumb screws (2). Lower seat pan.

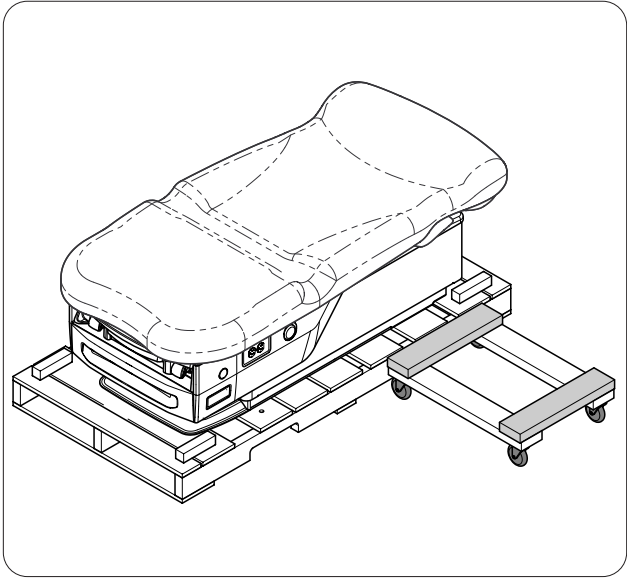
Step 6:
Place dolly next to head end of chair.

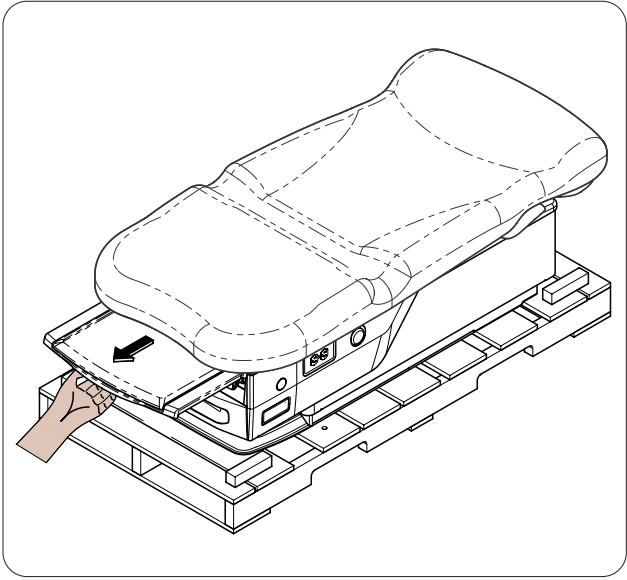
Step 7:
Slide footrest 9-10 inches out of chair.

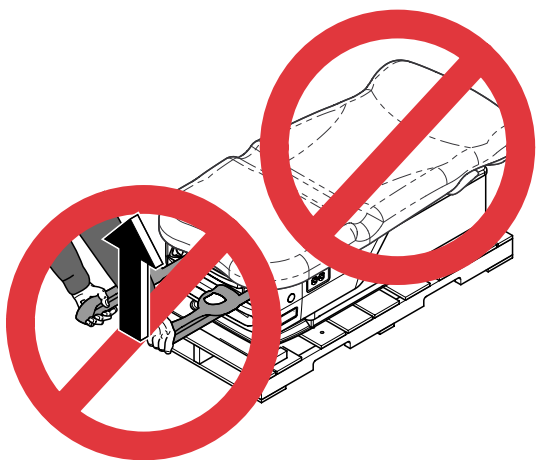
Equipment Alert
Do NOT lift chair by Upholstery
or Stirrups.

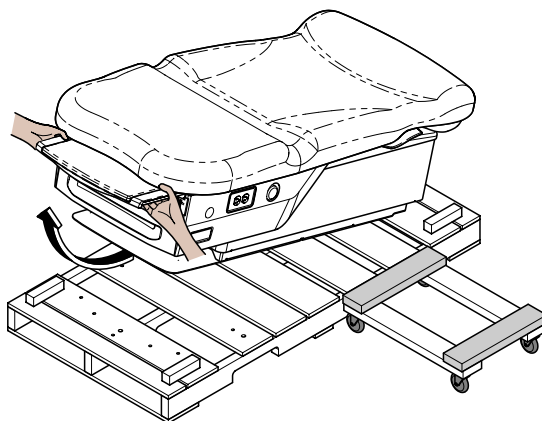
Step 8:
Using proper lifting technique, lift chair off front of skid.

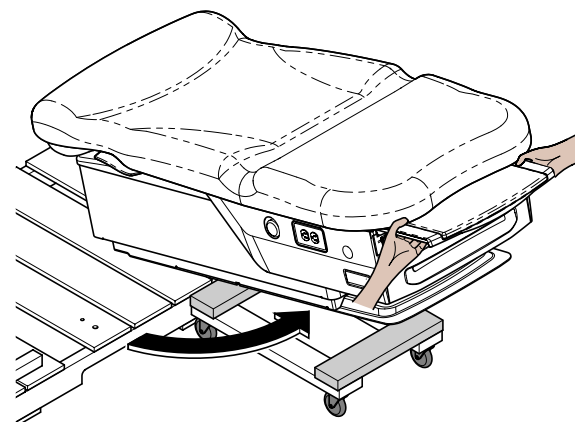
Step 9:
Rotate chair to dolly.

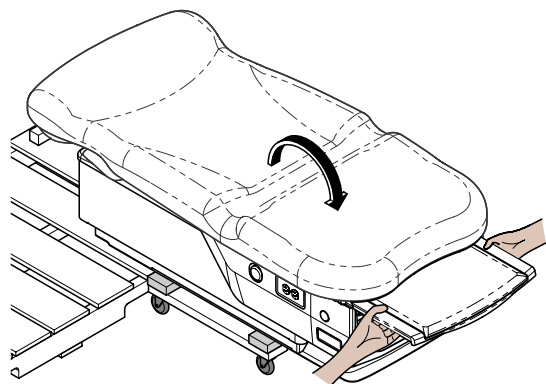Step 10:
Lower chair onto dolly.

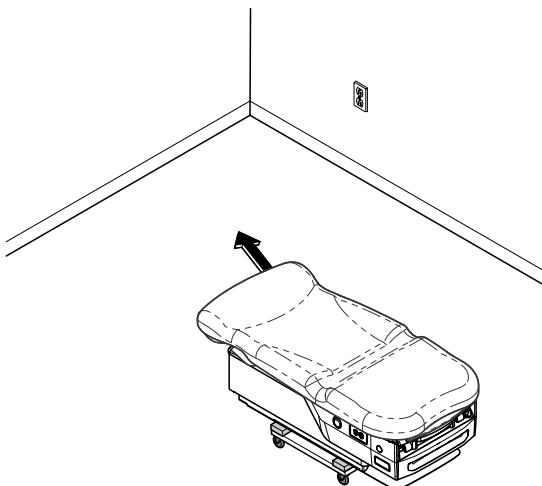
Step 11:
Move chair to desired location.

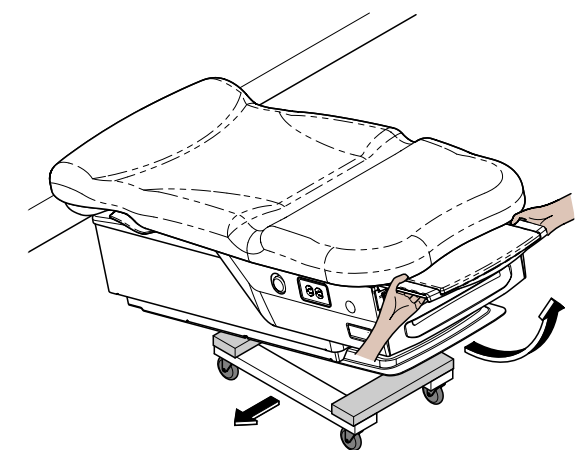
Step 12:
Remove chair from dolly.

Before connecting the Midmark device power cord into an electrical power source outlet, it is recommended to contact a local licensed electrician or a local authority having jurisdiction over the internal building power distribution to verify that the device when powered up will comply with all local, state and federal electrical codes. Use the electrical power rating marked on the device when determining the appropriate electrical branch circuit and outlet requirements. For 115V rated device models with a current rating at or above 12 Amps, a dedicated branch circuit may be required if the branch circuit protection is provided by a 15 Amp breaker.

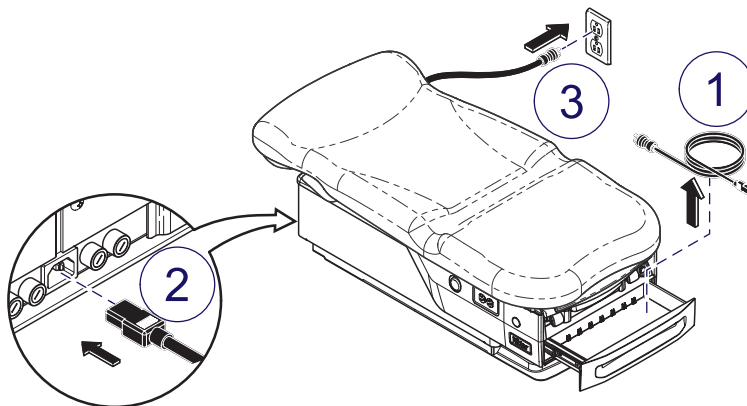
Step 13:

Remove power cord from drawer, connect to rear of chair and plug into wall or floor outlet.

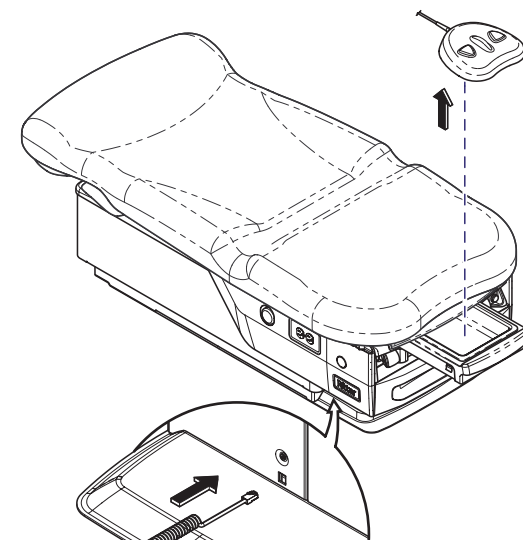
Step 14:

Remove foot control from treatment pan.
Connect foot control to front base of chair.

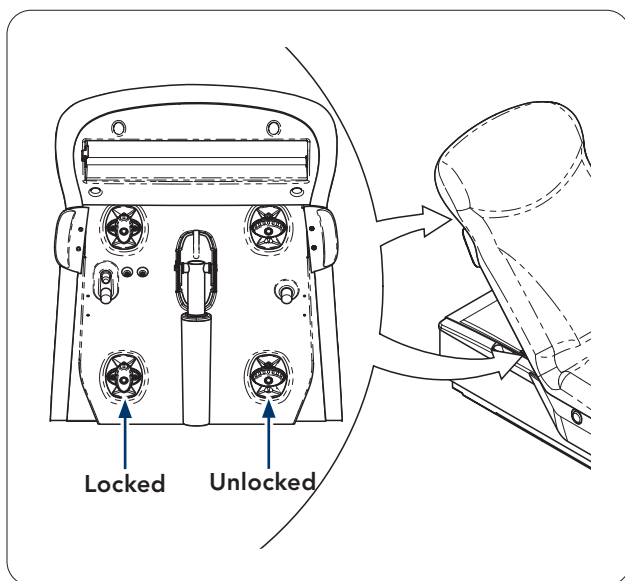Step 15:

Raise back and latch upholstery locking knobs (4).

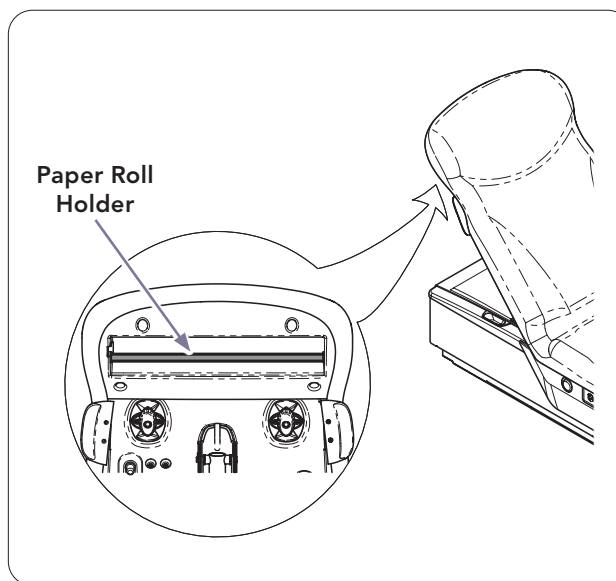
Step 16:

Ensure paper roll holder is in place.


Equipment Alert

The chair must be leveled or stable to ensure proper operation.

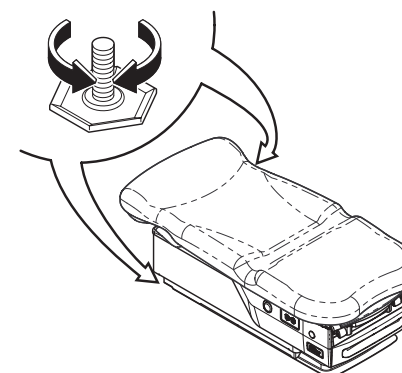
Step 17:


Adjust leveling screws (2) at head end of chair to level or balance chair.

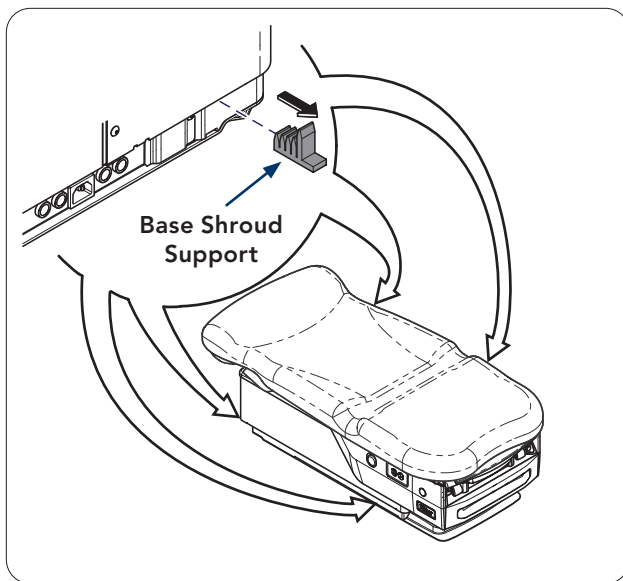Step 18:
Remove base shroud supports (4).

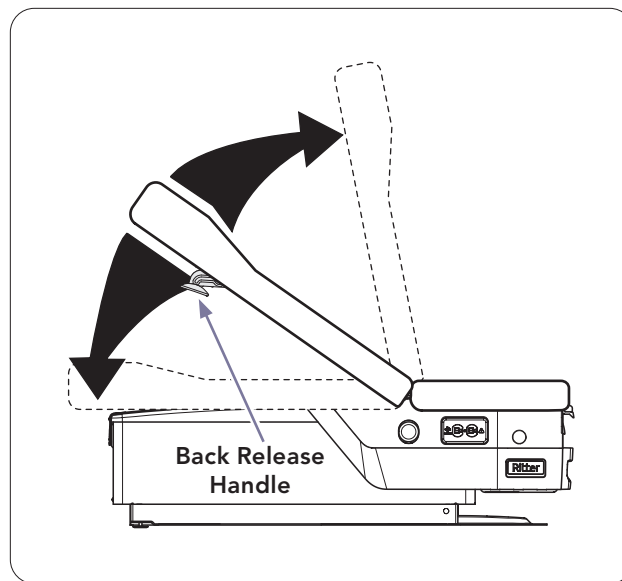
Step 19: (228 Only)
Verify proper manual Back operation.

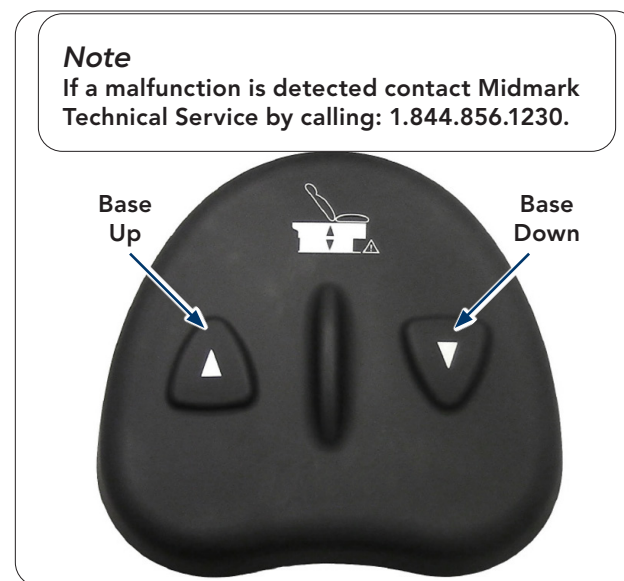
Note
If a malfunction is detected contact Midmark Technical Service by calling: 1.844.856.1230.

Step 20: (228 Only)
Verify proper operation of foot control. Press buttons for Base Up/Down.

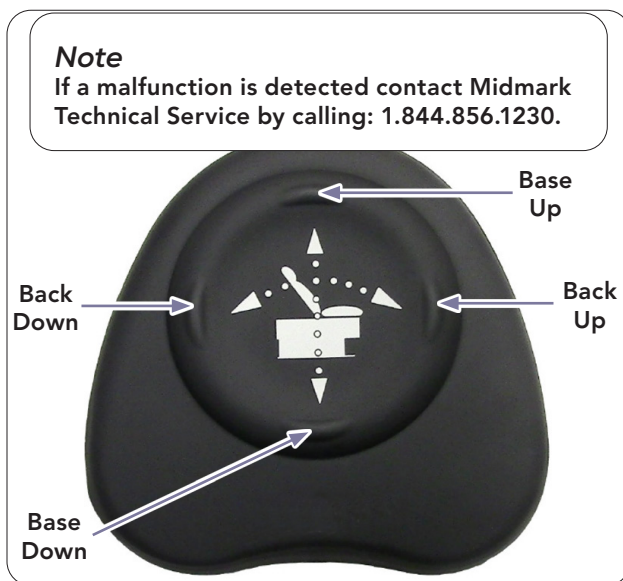Note
If a malfunction is detected contact Midmark Technical Service by calling: 1.844.856.1230.

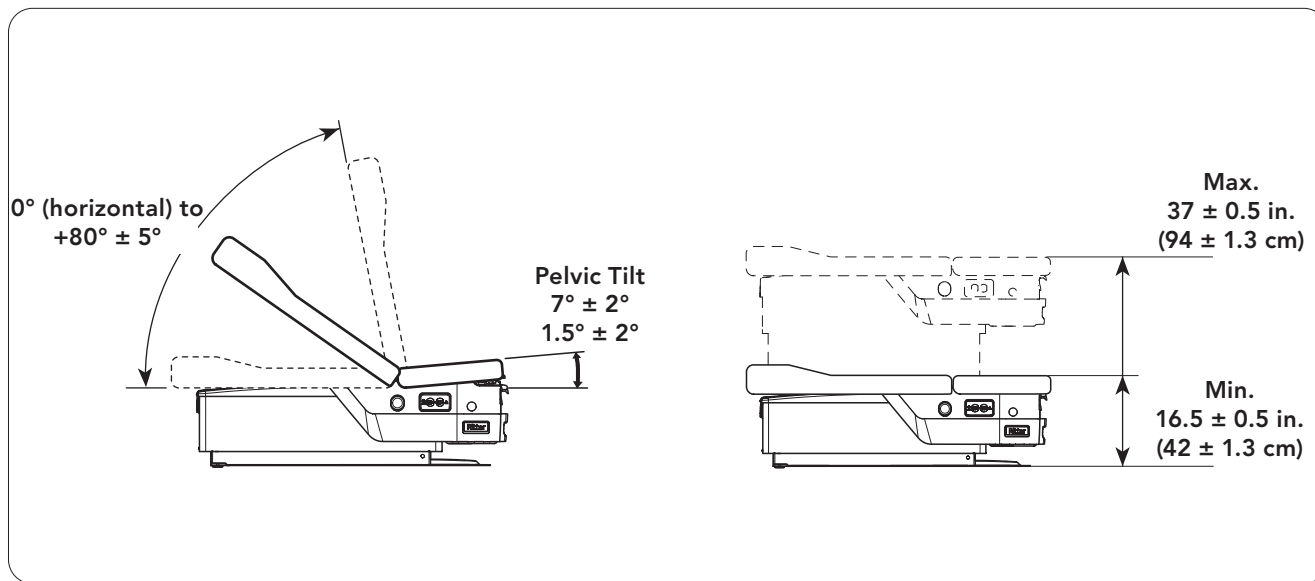Step 21: (229 Only)
Verify proper operation of foot control. Press buttons for Base Up/Down and Back Up/Down.

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