



Manual Examination Table

604



When you want a full-featured manual exam table, choose the Midmark 604

The Midmark 604 Manual Examination Table is designed to give you the best economy exam room option available with choice features. And this cost-effective table was designed to help enhance patient comfort. With the 604, the physician can easily position patients so they are comfortable and relaxed throughout the exam.

The exam table is where the practice of medicine takes place—the interaction between physician and patient

The new look of your exam room is intended to promote overall improvement of your office environment through the 604's stylized design, featuring a plush, one-piece upholstered top. The 604 not only looks great, but is designed to be simple to clean. In addition, the 604 is user-friendly with a pneumatic cylinder designed to assist in easily placing the counterbalanced top in an infinite number of positions. Release handles on both sides of the top section are designed to make positioning easy.



Ultra-Premium top option is more comfortable for your patients, giving your practice an exclusive feel and elevating your reputation.




**Exam Assistant®
Fold-Down Drawer System**
Provides an area to store basic pelvic and examination supplies at your fingertips. Its unique fold-down front can also turn into a work shelf for examination assistance. *(Patented)*



Flex Drawer Design
A flexible rear drawer design allows changing to either right-hand or left-hand access without the use of tools. It is designed to eliminate the potential for wall damage or the unkept look of pass-through drawer systems.



Premium Comfort System®
Unique elastomeric material stretched over a seat cut-out is designed to enhance patient comfort for a chair-like feel in an exam table. The upholstery is removable and available with optional heater. *(Patented)*

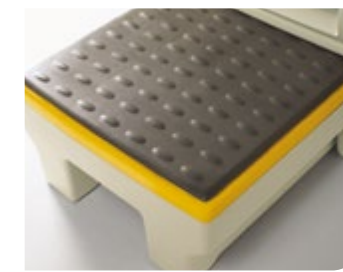

Large Foot Step
This rubberized footstep is designed to be easy to clean and slip resistant. The 40% larger footstep has bright yellow edging to assist the visually impaired in seeing the step edges clearly. *(Patented)*



Heavy Construction
The solid construction of the table is designed to support larger patients up to 500 pounds as well as provide many years of use.

Specifications

Height with Upholstery:
Premium/Ultra-Premium Top:
33.38" (84.8 cm)

Length with Footrest
Extended: 75.5" (191.8 cm)

Patient Weight Capacity:
500 lb (226 kg)

Electrical Receptacle:
(115 VAC, 5 AMPS Max.)
Duplex, hospital grade,
electrical receptacle
on patient's left side

Paper Roll Storage:
4 Rolls - 21" x 3.5"
(53.3 cm x 8.9 cm)

Paper Roll Holder:
Holds 21" x 3.5" Paper Roll
(53.3 cm x 8.9 cm)

Stirrups are standard

Electrical Ratings

115 VAC, 50/60HZ, 5.0 AMPS
(with receptacle)

115 VAC, 50/60HZ, 5.5 AMPS
(with drawer heater)

115 VAC, 50/60HZ, 6.5 AMPS
(with drawer heater and
heated top)

Safety Standards

UL 60601-1

CAN/CSA No. 601.1-M90

Options

**Premium/Ultra-Premium
Upholstery Top:**

28.5" W x 59.5" L
(72.4 cm W x 151.1 cm L)

Pelvic Tilt: 7°

Heated Drawer (45 watt)



Premium Solid Upholstery Colors

Blueberry-235

Cashmere Blue-239

Silver Sage-236

Terra Cotta-240

Fossil Grey-238

Tea Green-237

Black-312

Premium Pattern Upholstery Colors

Belagio-247

Capri-251

Siena-248

Treviso-252

Firenze-250

Tuscany-249

Ultra-Premium Solid Upholstery Colors

Ultra-Pacific Blue-241

Ultra-Sky-245

Ultra-Spruce-242

Ultra-Canyon-246

Ultra-Arctic-244

Ultra-Olivine-243

Ultra-Raven Wing-228



PC Work Shelf

Provides an adjustable height workspace for a tablet PC, portfolio or exam instruments and supplies.



Assist Handles

Designed to offer patients a stable platform for accessing the table and a sense of security during positioning.



Knee Crutches

A more ergonomic and comfortable option compared to standard stirrups. The ball joint design offers ultimate flexibility in patient positioning. Available in a non-articulating option.



IQvitals® Table Mount

Convenient IQvitals® swing arm table mount positions vitals at the side of the exam table to support efficient workflow and keeps cords secure and out of the way.



LED Exam Light

Easy access controls, high intensity light and an adjustable focal spot makes the Ritter 253 the perfect complement to any exam room. Available in table mount and mobile options.

CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and ISO 9001 Certified Company.
For more information or a demonstration, contact your Midmark dealer or call: 1-800-MIDMARK Fax: 1-800-365-2499
Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8392
or visit our website at midmark.com

© 2005 Midmark Corporation
Midmark Corporation, Dayton, OH.
Products subject to improvement changes without notice
Litho in U.S.A. 007-0562-00 Rev. H1 (11/16)

