
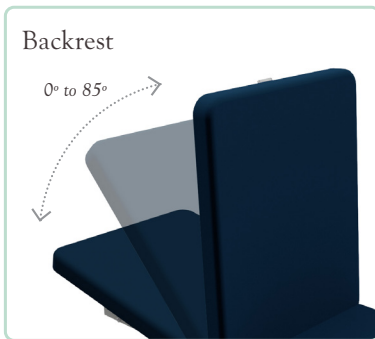


PT100/PT150 PHYSICAL THERAPY TABLES

Rock solid support for better results

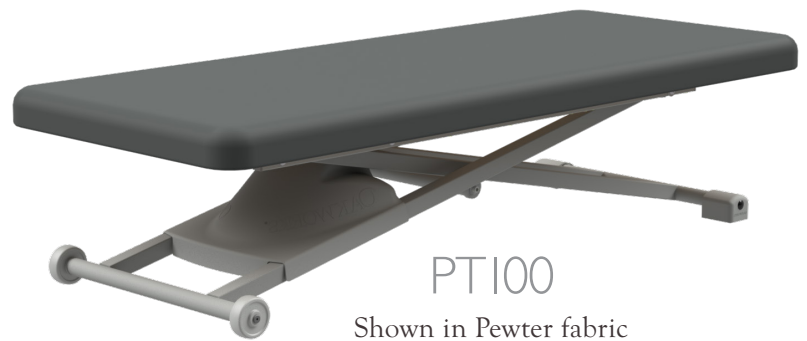

Features

- Super low 16" height
- Lift-Assist backrest 0° - 85° (PT150)
- Stable twin beam construction
- 550 lb weight capacity



Twin beam design lifts vertically, taking up much less space than cantilevered designs.

Visit www.oakworksmed.com for a quote!

PT150 2-SECTION TOP

- #72251-T* 29" W 1.75" Firm Response
- #73271-T* 31" W 1.75" Firm Response
- #71810-T* 29" W 2.5" Comfort
- #73402-T* 31" W 2.5" Comfort

PT100 RECTANGULAR TOP

- #68767-T* 29" W 1.75" Firm Response
- #73403-T* 27" W 1.75" Firm Response
- #68934-T* 29" W 2.5" Comfort
- #73404-T* 27" W 2.5" Comfort

FEATURES	
PT150 (29" & 31")	
2 Section Top	29" (74cm) or 31" (79cm) wide x 75" (191cm) long
Backrest	39" (99cm) long
Seat Section	36" (91cm) long
Motion 1	Electric 16" - 34" (40-86cm) (1.75" foam)
Motion 2	Manual Lift-Assist Backrest 0°-85°
PT100 (27" & 29")	
Rectangular Top	27" (68cm) or 29" (74cm) wide x 73" (185cm) long
Motion 1	Electric 17" - 34" (43-86cm) Height Range (1.75" foam)
Controls	
Controls	Foot control
Weight Capacity	550 lb (250kg)
Product Weight	151 lb (69kg)
Padding	1.75" (4cm) Firm Response or 2.5" (6cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	Accessory Mounting Plate (Backrest End)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Locking Twin Casters (one end of the table)	To easily move your table in the room
PT150 Side Rails	Drop Away (available for the 31" width)
Second Foot Control	Controls all powered movements (height)
Hand Control	Controls all powered movements (height)

ACCESSORIES <i>order any time.</i>	
Universal Paper Roll Holder	Mounts to installed accessory mounting plate
Face Hole Filler	Fits securely in cut out



Optional Locking Twin Casters


STABLE TWIN BEAM CONSTRUCTION

Stiffer than cantilever designs and with fewer moving parts, our unique design lifts vertically, taking up less room than a cantilevered design.




FIRM RESPONSE FOAM

Firm Response foam delivers a firm yet comfortable surface that maximizes ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design.

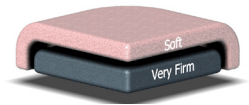

LIFT ASSIST BACKREST

Lift assist backrest features two powerful pneumatic cylinders that assist with adjusting the backrest up or down.

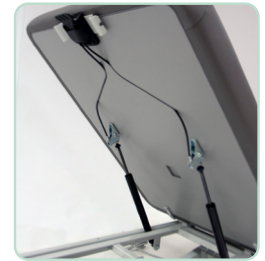

TERRATOUGH™ FABRIC



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

717.759.3146 www.oakworksmed.com email: medsales@oakworks.com