



PT250 PHYSICAL THERAPY TABLE

Rock solid support for better results



Shown in Pewter fabric

Features

- Super low 16" height
- Vertebral Artery Testing efficiently and safely
- Optimize Examination/Dry needling
- 550 lb lifting capacity



Large height range facilitates patient transfer & practitioner ergonomics.

PT250

#69193-T* 29" W 1.75" Firm Response #70234-T* 29" W 2.5" Comfort
 #68644-T* 31" W 1.75" Firm Response #70236-T* 31" W 2.5" Comfort

FEATURES	
2 Section Top	29" (74cm) or 31" (79cm) wide x 72" (183cm) long
Flat Section	61" (155cm) long
Head Section	11" (28cm) long
Motion 1	Electric Height Range 16" - 34" (40-86cm) (1.75" foam)
Motion 2	Manual position head section 40° to -90°
Controls	Foot control
Weight Capacity	550 lb (250kg)
Product Weight	151 lb (69kg)
Padding	1.75" (4cm) Firm Response or 2.5" (6cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	Accessory Mounting Plate (Foot End)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Locking Twin Casters (one end of the table)	To easily move your table in the room
Second Foot Control	Controls height
Hand Control	Controls height

ACCESSORIES <i>order any time.</i>	
Universal Paper Roll Holder	Mounts to installed accessory mounting plate
Face Hole Filler	Fits securely in cut out


STABLE TWIN BEAM CONSTRUCTION

Stiffer than cantilever designs and with fewer moving parts, our unique design lifts vertically, taking up less room than a cantilevered design.



FIRM RESPONSE FOAM

Firm Response foam delivers a firm yet comfortable surface that maximizes ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design.

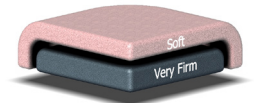

“The drop down head section feature enables patient safe and Physical Therapist efficient cervical examination and treatment.”

Dr. Robert DuVall,
 PT, DHSc, ATC, OCS, SCS, FAAOMPT CSCS



Optional Locking Twin Casters

TERRATOUGH™ FABRIC

Coal-T01	White-T26	Heron-T05	Earth-T07
Opal-T12	Stone-T13	Sapphire-T17	Clay-T18
Ocean-T20	Pewter-T29	Sage-T22	Blue Grass-T24
Orchid-T19	Espresso-T27	Sky Blue-T04	Ruby-T08

We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

717.759.3146 www.oakworksmed.com email: medsales@oakworks.com