



# PT400M PHYSICAL THERAPY TABLE

*Rock solid support for better results*

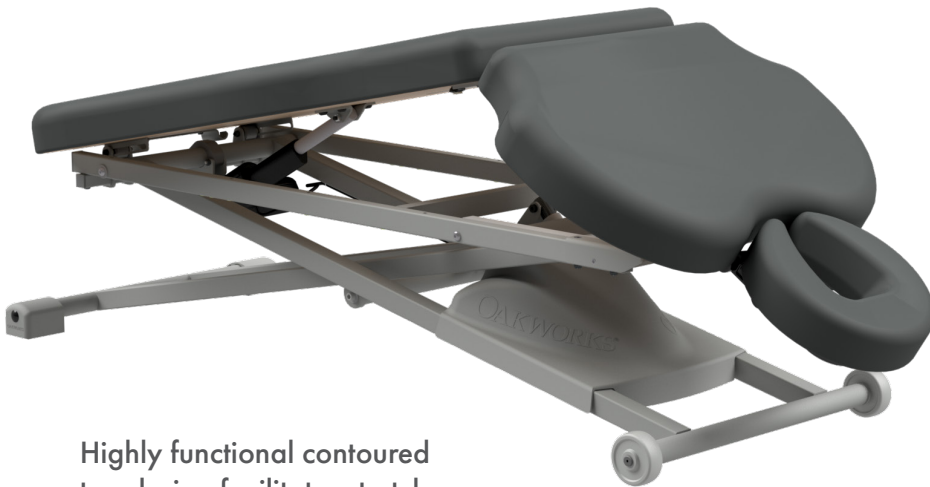

Highly functional contoured top design facilitates stretches and excellent ergonomics.

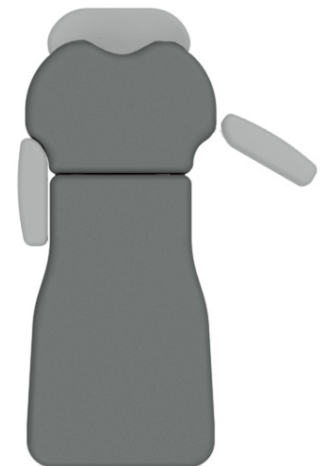

Optional Rotating Arm Rest and Removable Extension Pad as needed for support and patient access.

Powered Center Break  
0° to 45°



## Features

- Powered center break & height
- Contoured ergonomic design
- Advanced head support included
- Optional arm rests & table extenders



Twin beam design lifts vertically, taking up much less space than cantilevered designs.

# PT400M

#PKG73396-T\* 1.75" Firm Response

#PKG72175-T\* 2.5" Comfort Foam

FEATURES	
2 Section Contoured Top	29.5" (75cm) wide (at widest point) x 74" (188cm) long (including head rest)
Motion 1	Electric Height Range 17" - 34" (43-86cm) (1.75" foam)
Motion 2	Electric Center Break 0°-45°
Controls	Multi-function foot control
Weight Capacity	550 lb (250kg)
Product Weight	186 lb (84kg)
Padding	1.75" (4cm) Firm Response or 2.5" (6cm) Comfort Foam™
Fabric	Premium, PVC free TerraTouch™
Included Features	QuickLock™ Head Rest & Boiance™ pad, Accessory Mounting Plate (both ends)
Warranty	3 years parts, 2 years labor
Safety Listings	FDA listed, cETLus listed, CB Certificate

OPTIONS <i>must be ordered with table.</i>	
Locking Twin Casters (one end of the table)	To easily move your table in the room
Second Foot Control	Height and center break
Hand Control	Height and center break

ACCESSORIES <i>order any time.</i>	
Extension Pad	Use as table extender or lowered arm rest shelf (bracket required) 22.5" (57cm) wide x 9" (23cm) long
Rotating Arm Rests	Set of 2 4.5" (11cm) wide x 16" (41cm) long
Shelf Bracket	To use Extension Pad as a lowered arm rest shelf
Stabilization Strap	Stabilize your patient
Universal Paper Roll Holder	Mounts to installed accessory mounting plate
Face Hole Filler	Fits securely in cut out



Optional Locking Twin Casters


## STABLE TWIN BEAM CONSTRUCTION

Stiffer than cantilever designs and with fewer moving parts, our unique design lifts vertically, taking up less room than a cantilevered design.



## FIRM RESPONSE FOAM

Firm Response foam delivers a firm yet comfortable surface that maximizes ergonomic efficiency during examination and intervention procedures. Additionally, this surface minimizes Manipulation Dry Hits with its super high density, twin layer design.

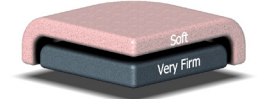

Removable Extension Pad & Shelf Bracket



Rotating Arm Rests removable accessory

## TERRATOUGH™ FABRIC



We strive to ensure that the colors displayed on our product imagery (online and in print) are as accurate as possible. However, due to the limitations of printers, scanners, computer monitors, and other equipment, the colors may not be accurate.

717.759.3146 [www.oakworksmed.com](http://www.oakworksmed.com) email: [medsales@oakworks.com](mailto:medsales@oakworks.com)