



A leader in infection prevention  
**since 1925**

**HIGH THERMAL DISINFECTION WASHER  
TIVA8-MD | 165 LITER CHAMBER**

Medical



Quality Products ♦ Leading Edge Technology ♦ Reliable Services ♦ Professional Excellence

# Welcome to the TIVA8-MD

## High Disinfection Thermal Washer

The ultimate automated solution for washing, rinsing, disinfecting, and drying your complex instruments prior to sterilization. Designed to simplify your work by eliminating manual washing and reducing your manpower. Ensuring the safety of your patients by providing consistent quality washing cycles at all times and ensuring the safety of your staff by eliminating risks of exposure to potentially harmful aerosols.



Ambulatory surgery centers



Hospitals CSSD




Physician's Office




Vet hospitals/clinics

## TIVA8-MD Washer Disinfecter Package

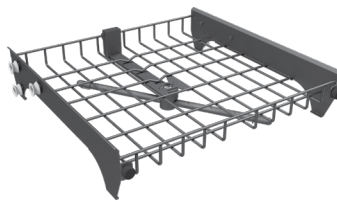

Upper Rack with spray arm

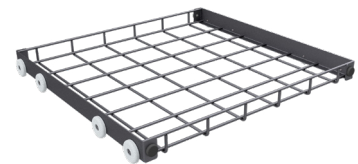

Lower Rack basic

### Features included:

- Double-layered tempered glass door
- Three water connections: Cold + Warm + Deionized water softener with protection filter
- Spray Arm Monitoring
- Multi Color Chamber LED Light
- Upper Rack with spray arm
- Lower Rack basic


### Comes in 2 different power configurations:

- 3-Phase 208-240V without Neutral, 60Hz  
Part# TIVA8-MD 3 MEDICAL
- Single Phase 208-240V without Neutral, 60Hz  
Part# TIVA8-MD 1 MEDICAL

**Medical Package comes standard with purchase of TIVA8 Medical Washers**



Upper Rack Support Grid 2 X DIN

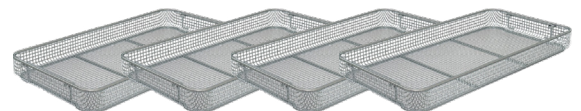

4 X Din Baskets (19" x 10" x 2")

# Improved Productivity and Efficiency



## Premium Performance

---

- Instruments are completely dried thanks to powerful electric blowers and heaters with a HEPA H14 filter
- Three water connections for flexible operation and for optimal washing needs (cold & hot with softeners, demineralized)
- Conductivity measurements to check rinse water quality before disinfection
- Up to 40 washing programs
- AISI 316L Stainless Steel



## Intuitive User Experience

---

- A full control and all-in-one operation from a single touch panel with color display
- Spray Arm monitoring for continuous washing validation
- Graphical report of washing phases including time and temperature
- Multi color chamber LED light (White, Green, Red) operates during loading, unloading and running the program
- Double-layered tempered glass door for visual monitoring



## Enhanced Traceability

---

- Comply with regulation by saving washing cycles data
- USB and RS232 connectivity
- Printer for data documentation (optional)




## Eco-Friendly

---

- 2 dosing pumps with flowmeter for precise consumption of acid and alkaline detergents
- Integrated water softener for reducing decontamination
- Built-in steam condenser for excess water and steam

# TIVA8-MD Technical Data

23.6 inch narrow design

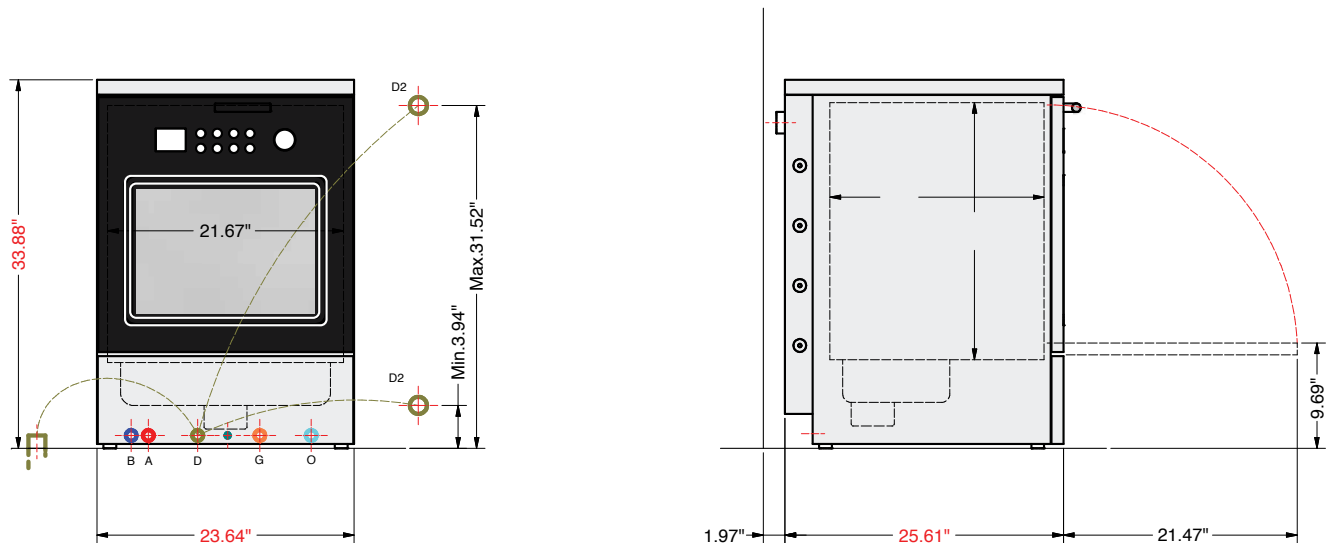

**A** Hot water connection-3/4" (tube 1/2" g)

**B** Cold water connection-3/4" (tube 1/2" g)

**D2** Drain dn 40 - hose connection 0.984" (at wall)

**G** Electrical connection

**O** Demineralized water connection-3/4" g

**E2** Drain dn 40 - hose connection 25mm (at wall)

## ALL YOUR INFECTION CONTROL NEEDS

Let us introduce you to our complementary family of sterile processing products to ensure a comprehensive process from start to end



**T-Clean TIVA Detergent:**  
Detergent for reprocessing thermostable and thermolabile instruments



**T-Clean TIVA Neutralizer:**  
Acidic neutralizing agent and detergent for the automated reprocessing of instruments and laboratory glassware



**T-Clean Prep-Foam:**  
Multi-Tiered Enzymatic Foam Spray. Helps To Prevent The Adhesion Of Bioburden To The Surfaces Of All Surgical Instruments And Scopes



T-Edge



Water Distiller



Biological Indicator Starter Kit

LTR330-0109/WAS/TIVA8MED/USA/V - 1/02.22

Medical

Distributed by:

**Tuttnauer**  
Innovation · Legacy · Partnership